

# 26" HD LED Medical Display

## FS-P2603D



The superb picture quality and wide viewing angles found in FS-P2603D offer the optimal visualization requirements demanded in today's OR. Medical professionals can view crystal clear images from a variety of locations in the room.




### Ultimate Compatibility

FS-P2603D accepts a large selection of video signal types, analog or digital. It also features rapid signal detection of imaging sources found in the OR. This means that FS-P2603D is compatible with nearly every surgical video integration application, old or new. It offers seamless deployment in surgical suites where many different video signals are present.

Signal inputs: DVI-D (2), 3G-SDI (2), SOG (2), VGA, C-Video, S-Video, Component Video.


Signal outputs: DVI-D (2), 3G-SDI (2), SOG (2), C-Video, S-Video, Component Video.

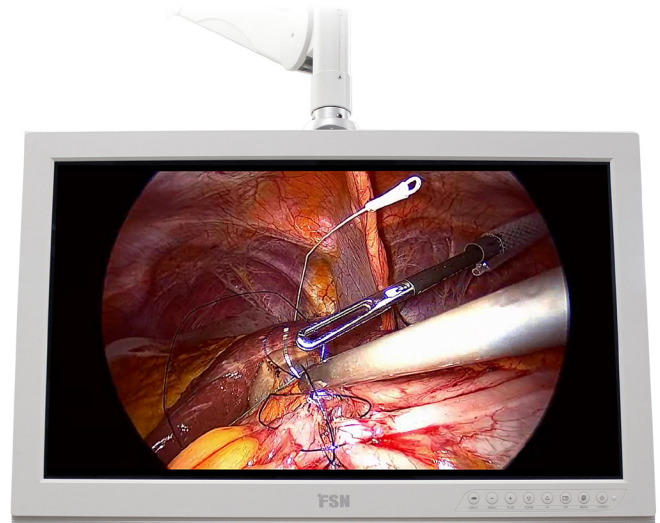

### Attention to Detail

FS-P2603D has the latest advanced features built in, including precise image controls and on-board DC out power for smaller component needs. Images can be zoomed, panned left or right, or freeze-frame. Two separate video signals can be shown as picture-in-picture or side-by-side.

FS-P2603D is constructed to medical-grade standards, and is ready out of the box for near-patient use in the sterile field. It can be easily mounted on spring arms, surgical booms, or mobile carts.

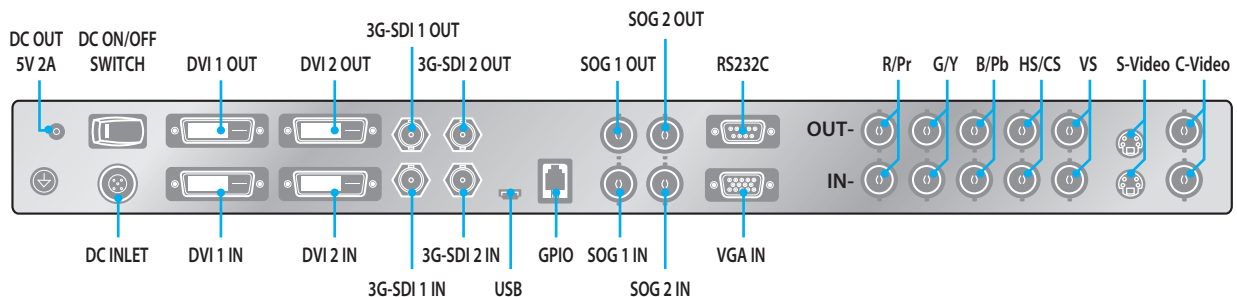

[www.fsnmed.com](http://www.fsnmed.com)

Foreseeson Custom Displays, Inc.  
2210 E. Winston Road  
Anaheim, CA 92806 USA  
Tel: 714-300-0540  
Fax: 714-300-0546

Latin America  
Aventura, FL 33180 USA  
Tel: 714-507-3855

Foreseeson Korea  
B-408, U-Space2  
670 Daewangpangyo-ro,  
Bundang-gu, Seongnam-si  
Gyeonggi-do, Republic of Korea  
Tel: +82-31-8017-0780

Foreseeson GmbH  
Industriestrasse 38a  
63150 Heusenstamm, Germany  
Tel: +49 6104 64398 0  
Fax: +49 6104 64398 11

Foreseeson UK Ltd.  
1 Wolsey Road  
East Molesay  
Surrey  
KT8 9EL  
United Kingdom  
Tel: +44 (0) 208 546 1047

Foreseeson (Shanghai)  
Medical Equipment Co., Ltd.  
Room 8E, No. 89 Building  
1122 North Qinzhou Road  
Xuhui, Shanghai 200233, China  
Tel: 86-21-6113-4188

# 26" HD LED Medical Display

## FS-P2603D






## Specifications

### General Features

Item	Description
Panel	26 inch TFT LCD (LED)
Resolution	1920 x 1080 pixel
Aspect Ratio	16 : 9
Active Area	576.0(H)mm x 324.0(V)mm
Pixel Pitch (mm)	0.30 x 0.30
Response Time	10 ms (typ. rise time)
Number of Colors	1.07 Billion
Brightness (typical)	450 cd/m <sup>2</sup>
Contrast Ratio (typical)	1400 : 1
Surface Treatment	Anti-glare
Viewing Angle (CR>10)	R/L 178°, U/D 178°
Input Signal	2 x DVI-D
	2 x SD/HD/3G-SDI (BNC)
	2 x SOG
	1 x VGA (D-sub)
	1 x C-Video (BNC)
Input Signal	1 x S-Video (DIN)
	1 x Component (RGBS, YPbPr) (5 x BNC)

Item	Description
Output Signal	2 x DVI-D
	2 x SD/HD/3G-SDI (BNC)
	2 x SOG
	1 x C-Video (BNC)
Output Signal	1 x S-Video (DIN)
	1 x Component (RGBS, YPbPr) (5 x BNC)
Power Supply	AC/DC Adaptor (AC 100~240V, DC 24V 6.25A)
Power Consumption	70W max
Accessories Included	User manual, AC-DC adaptor, AC power cord, DVI-D cable, BNC cable, Screw BH M4 x 14
Unit Dimension	668(W) x 421(H) x 88.8(D) mm 26.3(W) x 16.57(H) x 3.5(D) inch
Package Dimension	755.65(W) x 654.05(H) x 234.95(D) mm 29.75(W) x 25.75(H) x 9.25(D) inch
IP Rating	IP32 - overall, IP44 front side only
IK Rating	IK05 with 3H AR glass
Weight	8.63 kg, 19.03 lbs. (monitor with cover)
	13.2 kg, 29.1 lbs. (shipping package)
Certifications	ANSI/AAMI ES60601-1 (UL60601-1), CAN/CSA C22.2 No. 60601-1, FCC Class B, EN60601-1, EN60601-1-2, CE, MDR(EU) 2017/745 Class I Medical Device

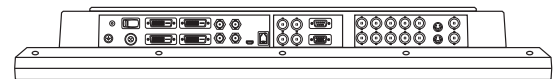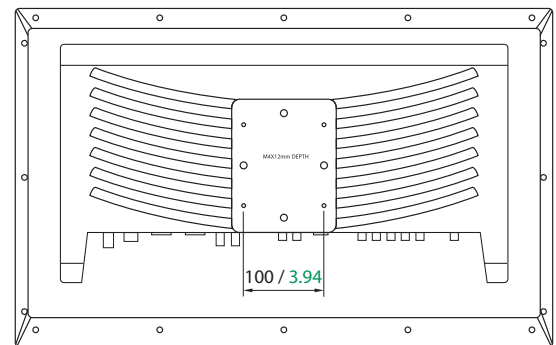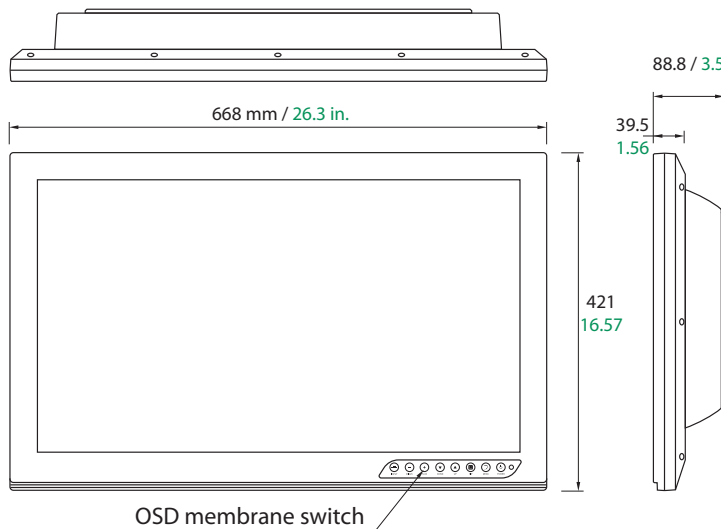

FSN Medical Technologies products meet EU RoHS Directive 2011/65/EU+2015/863/EU and Registration, Evaluation, Authorisation and Restriction of Chemicals (REACH) obligations. FSN manages the information flow as required by regulations.

Specifications are subject to change with or without notice. Doc. # FSN1948 Rev. 11/21

[www.fsnmed.com](http://www.fsnmed.com)

#### Foreseeson Custom Displays, Inc.

2210 E. Winston Road  
Anaheim, CA 92806 USA  
Tel: 714-300-0540  
Fax: 714-300-0546

Latin America  
Aventura, FL 33180 USA  
Tel: 714-507-3855

#### Foreseeson Korea

B-408, U-Space2  
670 Daewangpangyo-ro,  
Bundang-gu, Seongnam-si  
Gyeonggi-do, Republic of Korea  
Tel: +82-31-8017-0780

#### Foreseeson GmbH

Industriestrasse 38a  
63150 Heusenstamm, Germany  
Tel: +49 6104 64398 0  
Fax: +49 6104 64398 11

#### Foreseeson UK Ltd.

1 Wolsey Road  
East Molesey  
Surrey  
KT8 9EL  
United Kingdom  
Tel: +44 (0) 208 546 1047

#### Foreseeson (Shanghai)

Medical Equipment Co., Ltd.  
Room 8E, No. 89 Building  
1122 North Qinzhou Road  
Xuhui, Shanghai 200233, China  
Tel: 86-21-6113-4188