


26" HD LED Medical Display

FS-P2604D



Designed for a Better Image

This 26 inch medical grade display from FSN has the latest technology needed to reproduce super-sharp HD images. The front filter bonding process provides higher brightness, higher contrast ratio, minimized glare, resistance to internal condensation, and increased strength. Flush-front glass

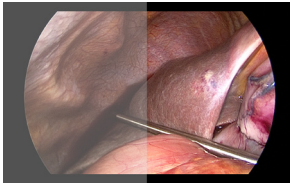


FSN's bonding process (right) improves picture quality.

with touch interface control is easy to clean and allows for fast procedure turnaround. FS-P2604D has an ingress protection (IP) rating of 33 to safeguard the display monitor's internal components and protect against the intrusion of solids or liquids.

Broad Compatibility

Thanks to rapid signal detection and robust mode tables, FS-P2604D automatically recognizes other highly specialized surgical and diagnostic equipment found in surgical suites, operating rooms, emergency rooms, and procedural facilities.

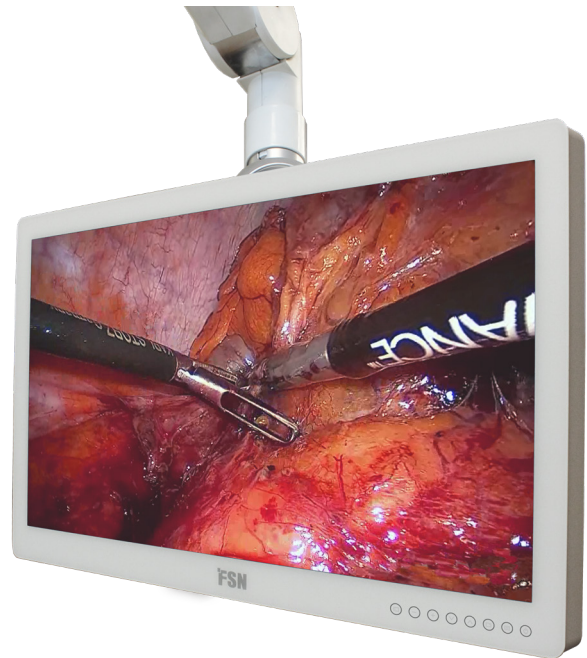


Video connections are available for DVI-D, SDI (SD/HD/3G-SDI), VGA, C-Video, S-Video, and Component. This ensures compatibility with popular signal types found in surgical applications.



User Friendly

FS-P2604D features an intuitive, full-featured on screen menu for image and layout settings. Multi-modality options let doctors view two images at the same time. The rounded monitor housing corners help protect OR personnel from accidental bumping injuries. Its lightweight design is ideal for boom or spring-arm mounting, and the convenient back cable cover helps organize and manage multiple cords.

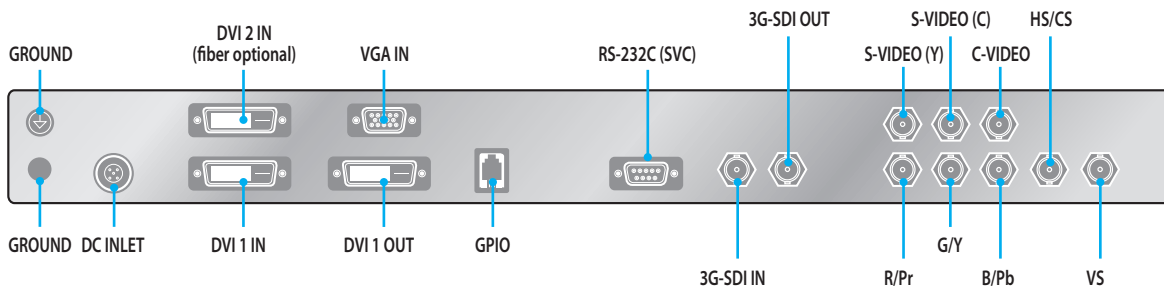

www.fsnmed.com

Foreseeson Custom Displays, Inc.

2210 E. Winston Road
Anaheim, CA 92806 USA
Tel: 714-300-0540
Fax: 714-300-0546

Latin America
Aventura, FL 33180 USA
Tel: 714-507-3855

Foreseeson Korea

B-408, U-Space2
670 Daewangpangyo-ro,
Bundang-gu, Seongnam-si
Gyeonggi-do, Republic of Korea
Tel: +82-31-8017-0780

Foreseeson GmbH

Industriestrasse 38a
63150 Heusenstamm, Germany
Tel: +49 6104 64398 0
Fax: +49 6104 64398 11

Foreseeson UK Ltd.

1 Wolsey Road
East Molesay
Surrey
KT8 9EL
United Kingdom
Tel: +44 (0) 208 546 1047

Foreseeson (Shanghai)

Medical Equipment Co., Ltd.
Room 8E, No. 89 Building
1122 North Qinzhou Road
Xuhui, Shanghai 200233, China
Tel: 86-21-6113-4188

26" HD LED Medical Display

FS-P2604D







Specifications

General Features

Item	Description
Panel	26 inch TFT LCD (LED)
Resolution	1920 x 1080 pixel
Aspect Ratio	16 : 9
Active Area	576.0(H)mm x 324.0(V)mm
Pixel Pitch (mm)	0.30 x 0.30
Response Time (typical)	10 ms (rise time)
Number of Colors	1.07 Billion
Brightness (typical)	500 cd/m ²
Contrast Ratio (typical)	1500 : 1
Surface Treatment	Anti-glare
Viewing Angle (CR>10)	R/L 178°, U/D 178°
Input Signal	1 x DVI-D 1 1 x DVI-D 2 (fiber optional) 1 x VGA (D-sub) 1 x SD/HD/3G-SDI (BNC) 1 x C-Video (BNC) 1 x S-Video (Y/C) (2 x BNC) 1 x Component (RGBS, YPbPr) (5 x BNC)

Item	Description
Output Signal	1 x DVI-D 1 x SD/HD/3G-SDI (BNC)
Power Supply	AC/DC Adaptor (AC 100~240V, DC 24V 6.6A)
Power Consumption	65W max
Accessories Included	User manual, AC-DC adaptor, AC power cord, DVI-D cable, BNC cable, Screw BH M4 x 10
Unit Dimension	640(W) x 402(H) x 78(D) mm 25.20(W) x 15.83(H) x 3.07(D) inch
Package Dimension	755.65(W) x 654.05(H) x 234.95(D) mm 29.75(W) x 25.75(H) x 9.25(D) inch
IP Rating	IP33 - overall, IP44 - front side only
IK Rating	IK07 with 5H AR glass
Weight	8.44 kg, 18.61 lbs. (monitor with cover) 12.85 kg, 28.33 lbs. (shipping package)
Certifications	ANSI/AAMI ES60601-1 (UL60601-1), CAN/CSA C22.2 No. 60601-1, FCC Class B, EN60601-1, EN60601-1-2, CE, MDR(EU) 2017/745 Class I Medical Device

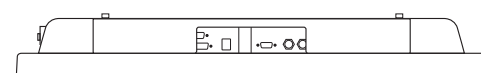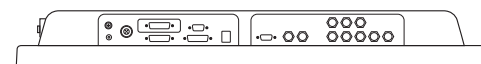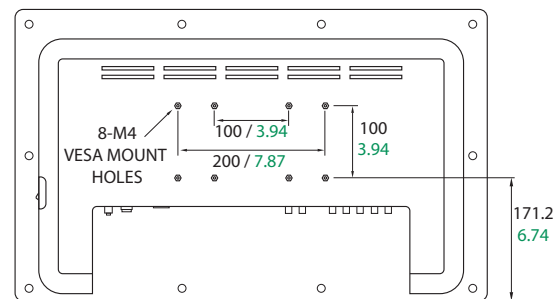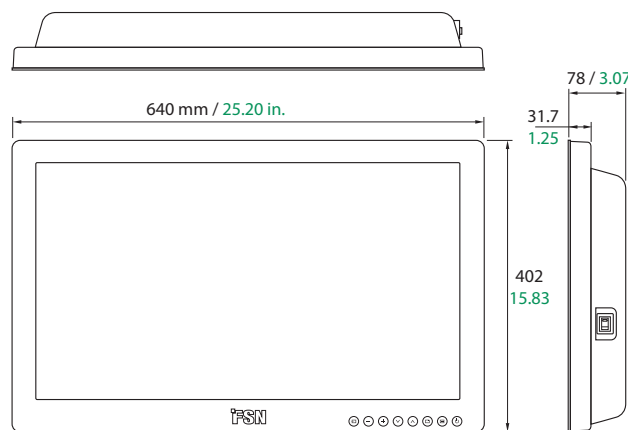

FSN Medical Technologies products meet EU RoHS Directive 2011/65/EU+2015/863/EU and Registration, Evaluation, Authorisation and Restriction of Chemicals (REACH) obligations. FSN manages the information flow as required by regulations.

www.fsnmed.com

Specifications are subject to change with or without notice. Doc. # FSN1987 Rev. 11/21

Foreseeson Custom Displays, Inc.

2210 E. Winston Road
Anaheim, CA 92806 USA
Tel: 714-300-0540
Fax: 714-300-0546

Latin America
Aventura, FL 33180 USA
Tel: 714-507-3855

Foreseeson Korea

B-408, U-Space2
670 Daewangpangyo-ro,
Bundang-gu, Seongnam-si
Gyeonggi-do, Republic of Korea
Tel: +82-31-8017-0780

Foreseeson GmbH

Industriestrasse 38a
63150 Heusenstamm, Germany
Tel: +49 6104 64398 0
Fax: +49 6104 64398 11

Foreseeson UK Ltd.

1 Wolsey Road
East Molesay
Surrey
KT8 9EL
United Kingdom
Tel: +44 (0) 208 546 1047

Foreseeson (Shanghai) Medical Equipment Co., Ltd.

Room 8E, No. 89 Building
1122 North Qinzhou Road
Xuhui, Shanghai 200233, China
Tel: 86-21-6113-4188